




IC-VIS Colorimeter

02/2025, V1.2

Safety manual

English 3

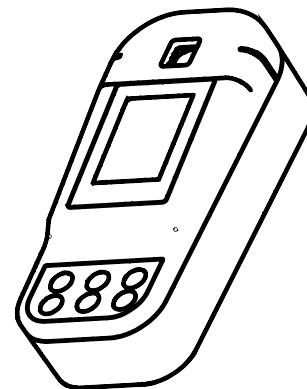

CONTENT

1	General information	5
1.1	Technical description.....	5
1.2	Safety.....	5
1.2.1	Symbols.....	5
1.2.2	Safety instructions	5
1.3	Compliance and certification	6
1.4	CE-marking	6
1.5	ICES-003 Compliance (Class A)	7
1.6	FCC Part 15 Compliance (Class A)....	8
1.7	Product label	9
2	Specifications.....	9
3	Installation	11
3.1	Product overview.....	11
3.2	Scope of delivery.....	12
3.3	Battery installation	13
3.4	Install sample holder.....	15
4	UI and navigation	16
4.1	Keypad	16
4.2	Display.....	17
5	Start-up	17
5.1	Turn device on	17
5.2	Configuration	18
5.2.1	General settings	18
5.2.2	Date-time	18
6	Standard operation.....	19

IC-VIS **SAFETY MANUAL** v1.2
Spec-Imaging 2025

6.1	Select method	19	7.2	Storage	26
6.1.1	Via menu	19	7.3	Battery replacement	26
6.1.2	Favourites	19	8	Troubleshooting	26
6.2	Add new methods	20	9	Service	27
6.2.1	Method from method library ..	20	9.1	Replacement parts and accessories	27
6.2.2	Method from calibration	21	9.2	Disposal	27
6.3	Perform measurement	21	9.3	Warranty, liability and complaints	27
6.4	Procedure for correct values	22		WARRANTY SUMMARY FOR SPEC-IMAGING	
6.5	IC – Instant Clarity turbidity compensation	23		PRODUCTS	27
6.6	Turbidity measurement	24	9.4	Contact	29
6.6.1	Turbidity calibration	24	9.5	Version history	29
7	Maintenance	25			
7.1	Instrument cleaning	25			

1 GENERAL INFORMATION

1.1 TECHNICAL DESCRIPTION

The IC-VIS is a portable colorimeter that enables attenuation and turbidity measurements at the wavelengths of 420 nm, 520 nm, 560 nm and 610 nm generated by a LED light source in combination with interference filters. The attenuation of light through the liquid sample is converted to a concentration of a certain chemical component, depending on the reagent used. Turbidity (as Formazine turbidity unit, FTU) is measured simultaneously with all measurements.

1.2 SAFETY

1.2.1 Symbols

Symbol	Description
	This symbol indicates that there is a risk for fire present.
	This symbol indicates that there is a risk for chemical exposure and that only authorized personnel should perform tasks in relation to equipment maintenance and handling of chemicals.
	This symbol indicates explosion hazard.

1.2.2 Safety instructions

Below are the 4 used hazard warnings.

 DANGER

Indicates an imminent or potentially hazardous situation which, if not avoided, will result in death or serious injury.

 WARNING

Indicates an imminent or potentially hazardous situation which, if not avoided, could result in death or serious injury.

 CAUTION

Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury.

NOTICE

Indicates a situation which, if not avoided, may result in malfunction or damage of the instrument. Information that clarifies facts

in the text and that requires special attention.

1.3 COMPLIANCE AND CERTIFICATION

This section outlines the certifications and regulatory compliances applicable to the IC-VIS. Adherence to these standards ensures that the product meets safety, performance, and environmental requirements as mandated by international and regional authorities.

1.4 CE-MARKING



IC-VIS **SAFETY MANUAL** v1.2
Spec-Imaging 2025

The CE marking declares that the product complies with the following European Community harmonisation legislation:

Electromagnetic Compatibility (EMC) Directive 2014/30/EU.

Restriction of Hazardous Substances (RoHS) Directive 2011/65/EU and 2015/863.

Directive 2012/19/EU on waste electrical and electronic equipment (WEEE).

This product fulfills the requirements for RF emission according to the Class A limit of EN 61326-1. Notice! This is a Class A product. In a domestic environment this product may cause RF interference, in which case the user may be required to take adequate measures.

You can obtain a copy of the original declaration of conformity from Spec-Imaging AB. Contact information can be found in Section 9.4.

1.5 ICES-003 COMPLIANCE (CLASS A)

Canadian Radio Interference-Causing Equipment Regulation, ICES-003, Class A:

Supporting test records reside with the manufacturer. This Class A digital apparatus meets all requirements of the Canadian Interference-Causing Equipment Regulations.

Cet appareil numérique de classe A répond à toutes les exigences de la réglementation canadienne sur les équipements provoquant des interférences.

1.6 FCC PART 15 COMPLIANCE (CLASS A)

Supporting test records reside with the manufacturer. The device complies with Part 15 of the FCC Rules. Operation is subject to the following conditions:


1. The equipment may not cause harmful interference.
2. The equipment must accept any interference received, including interference that may cause undesired operation.

Changes or modifications to this equipment not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This equipment has been tested and found to comply with the limits for a Class A digital

device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference, in which case the user will be required to correct the interference at their expense. The following techniques can be used to reduce interference problems:

1. Disconnect the equipment from its power source to verify that it is or is not the source of the interference.
2. Move the equipment away from the device receiving the interference.
3. Reposition the receiving antenna for the device receiving the interference.
4. Try combinations of the above.

1.7 PRODUCT LABEL

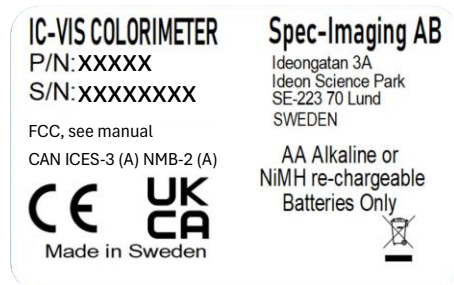

2 SPECIFICATIONS

The IC-VIS specifications are listed below

Specifications	Details
Battery	4 AA, around 2 months @ 10 measurements/day (5 days/week)
Certifications	CE, UKCA, FCCm CAN ICES-3 (A) NMB-2 (A)
Content	Colorimeter, unit lid, cuvette holders, sample cover, blank cuvettes, batteries, tool, start-guide
Cuvette compatibility	10 mm square, 13 and 16 mm round

Data logging/transfer	100 latest measurements/QR-code for data transfer
Detector	CMOS
Dimensions	215 mm x 96 mm x 60 mm
Display	LCD-IPS, backlit
Language	English, French, Swedish
Light source	LED
Measurement type	Concentration, absorbance/optical depth/ transmission, turbidity
Methods	Preprogrammed and programmable
Operating temperature	0-40 °C (max relative humidity 90% non-

	condensing; max 2000m)
Photometric accuracy/ range	± 0.03 abs / 0-3 ^a abs
Photometric repeatability (std)	± 0.001 ^b abs
Power requirement	4 AA LR6 alkaline batteries or HR6 NiMH AA rechargeable batteries
Protection class	Class III
Turbidity range	0-600 FTU
Turbidity compensation range	0-400 ^c FTU/ particles >5 μ m
Turbidity accuracy	± 10 FTU ^d
Storage temperature	-25-48 °C (max relative humidity 90% non-condensing)
Warranty	1 year

Wavelength selection	Method dependent - automatic
Wavelengths/Bandwidth	420 nm/10nm, 520 nm/10nm, 560 nm/10nm, 610 nm/10 nm
Weight	690 g

a – Up to 3 abs for 520, 560 and 610 nm. Up to 2 abs for 420 nm.

b- at 520 nm

c – for measurements done at 520, 560 and 610 nm. For measurements done at 420 nm the range is 0-200 FTU.

d - for turbidity measurements at 420 nm.

3 INSTALLATION

3.1 PRODUCT OVERVIEW

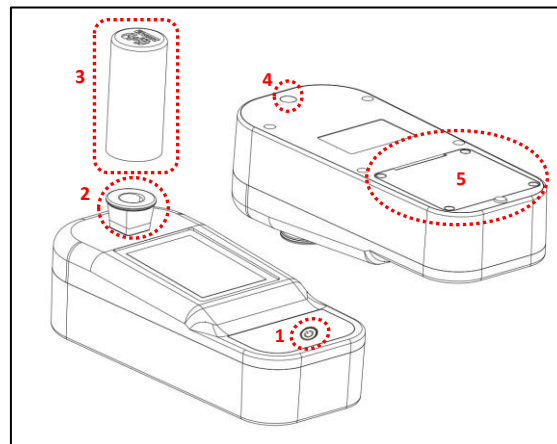


Figure 1. IC-VIS overview.

The IC-VIS is shown in Figure 1 with indicated details, listed below.

1	Power button	4	QR-code camera
2	Sample compartment	5	Battery compartment
3	Universal lid		

WARNING

It is under no circumstances allowed to open the device enclosure beyond battery replacement, as operator may be exposed to a strong light source or moving parts. Opening the enclosure is considered as improper use of the device and voids the warranty of the device.

3.2 SCOPE OF DELIVERY

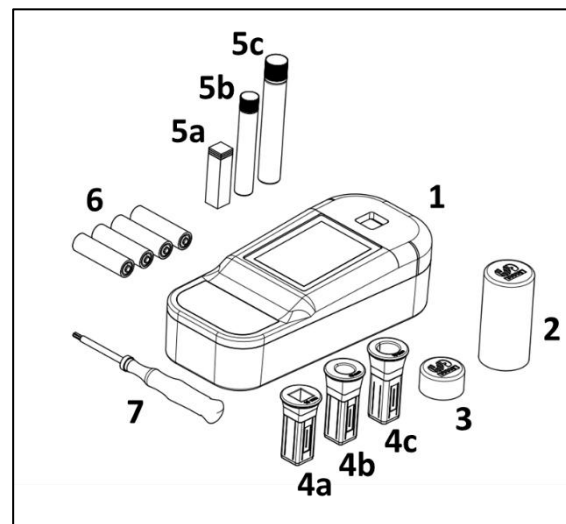


Figure 2. Overview of all parts delivered.

The components included in the package are shown in Figure 2 and listed below. Please

contact Spec-Imaging if any components are missing upon delivery.

1. IC-VIS Colorimeter
2. Sample Cover
3. Unit Cover
4. Cuvette holders
 - a. 1x 10mm square
 - b. 1x 13mm round
 - c. 1x 16mm round
5. Blank cuvettes
 - a. 2x 10 mm
 - b. 2x 13 mm
 - c. 2x 16 mm
6. 4 AA LR6 Alkaline batteries
7. Torx T10 screwdriver
 - Cleaning tissue for cuvettes
 - Quick start guide

3.3 BATTERY INSTALLATION

 WARNING	
	Fire hazard: Use only battery types specified to avoid risk of fire.

 WARNING	
	Explosion hazard: Ensure that batteries are of specified type and installed in the correct orientation. Mixing old and new batteries are to be avoided.

 CAUTION	
Use LR6 Alkaline AA or HR6 NiMH AA rechargeable batteries exclusively. Never use FR6 Lithium AA batteries.	

⚠ CAUTION

Make sure to indicate the battery type in the IC-VIS settings for optimal battery lifetime and device performance.

The Spec-Imaging product is powered by four AA batteries, which can be replaced by the user as needed. Make sure that all batteries are oriented in the correct direction as indicated in the battery compartment. See Figure 3.

For optimal battery lifetime and accurate device performance, make sure to set the correct battery type in the device settings. Main Menu → Settings → Battery.

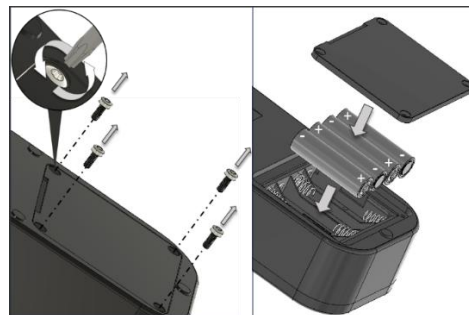


Figure 3. Battery installation.

The Spec-Imaging product uses a 3.0 V CR2032 lithium battery as the power supply for its internal real-time clock (RTC).

⚠ CAUTION

Do not replace the RTC battery yourself. Contact Spec-Imaging support if there is need for a battery change.

3.4 INSTALL SAMPLE HOLDER

The device has 3 types of sample holders (for 10 mm square, 13 and 16 mm round cuvettes). The cuvette should be installed with the cuvette size information in the correct direction (e.g. 13 mm), see Figure 4.

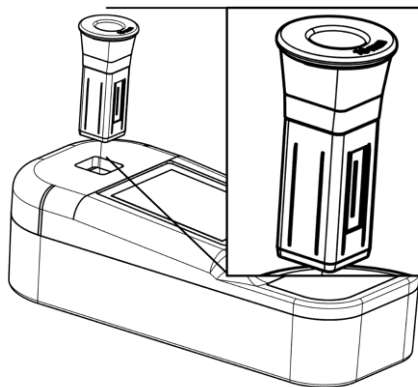


Figure 4. Cuvette holder orientation.

The sample holder needs to be inspected with regular intervals to make sure it is in good condition for accurate measurement procedure.

NOTICE

The correct orientation of the cuvette holders are crucial for correct operation of the device.

To swap cuvette holder grab hold of the cuvette holder top and pull.

4 UI AND NAVIGATION

4.1 KEYPAD

The keypad consists of 6 keys shown in Figure 5 and described below.



Figure 5. Overview of the button pad.



Power on / power off*.



Go to MENU. If in Menu, press for returning to main screen. Read QR-code*.



For BLANK measurement. If in menu, press for BACK.



For MEASUREMENT of sample after blank is performed. If in menu, press to SELECT.



UP if in list mode. Change unit if in main screen.



DOWN if in list mode. Enter FAVORITES LIST in Main Screen. Swap sub-metod* if in main screen.

*longpress

4.2 DISPLAY

The main screen is shown in Figure 6 and detailed description is outlined below.

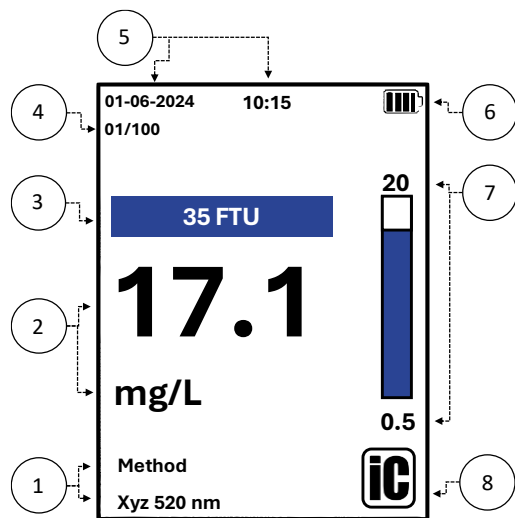


Figure 6. Overview of the main screen.

- 1- Method ID and parameter
- 2- Measurement value and unit
- 3- Turbidity reading
- 4- Measurement number
- 5- Date and time
- 6- Battery indicator
- 7- Method range
- 8- iC mode status

5 START-UP

5.1 TURN DEVICE ON

When turning on the unit, press the power button  and wait until main screen appears. Start screen is shown in Figure 6 with detailed description of indicators in section 4.2.

When turning off the device, put on the unit cover and do not store the device with a cuvette in the cuvette holder.

5.2 CONFIGURATION

The basic configuration of the device can be done by visiting the IC-VIS configuration page www.spec-imaging.com/ic-vis-configurator.



5.2.1 General settings

The general settings of the device can be set on the above specified IC-VIS configuration page.

Language: ENG

Sleep time: 5/10/20/30/45/60/90/120 min

*Number of measurement points in QR-code exporting:
1-20 points*

Color theme

Chose the preferred settings and scan the displayed QR-code by long-pressing the menu button . Confirmation of the new settings will be displayed as a message on the screen.

5.2.2 Date-time

Set the date and time of the device by visiting the IC-VIS configuration page and select **Set**

date/time. Longpress the menu button  to scan the QR-code. Make sure the correct date and time is set on the top of the main screen.

6 STANDARD OPERATION

The device is delivered with a list of 4 basic methods and can be further configured to contain up to 46 additional methods of any method provider.

6.1 SELECT METHOD

Measurements of a sample need to be done using the correct method which can be chosen in the following ways.

6.1.1 Via menu

A method can be selected by entering the list of methods in main menu → select method.

6.1.2 Favourites

Favourite methods list can be accessed by pressing DOWN , see Figure 7. Favourites can then be stepped through using the arrows  and , and selected with SELECT .

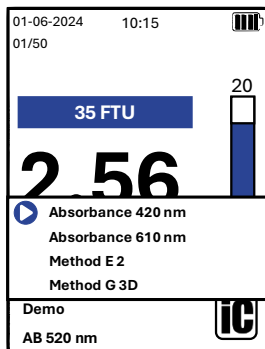


Figure 7. Favourites list shown when pressing DOWN when in main screen.

6.1.2.1 Add to favourites

Methods can be added to the favourites list by viewing a method in the methods list and long-pressing SELECT  until a star appears to the right of the method. The same procedure can be repeated to remove a method from the favourites list.

6.2 ADD NEW METHODS

New methods can be added to a device by importing them from the IC-VIS pages.

6.2.1 Method from method library

There are pre-calibrated methods in the **Methods library** which can be found on the IC-VIS Configurator page www.spec-imaging.com/method-library.



6.2.2 Method from calibration

It is also possible to do a method calibration manually and then adding it as a method in your device, through the IC-VIS Configurator page www.spec-imaging.com/ic-vis-configurator.



6.3 PERFORM MEASUREMENT

A measurement is performed by following the instructions on the main screen.

Insert blank – insert a blank sample, cover the sample with the lid and press BLANK .

NOTICE

Make sure to use the correct type of blank for your method. Sometimes a method blank is needed for accuracy.

Insert sample – insert the sample (prepared according to method specific description), cover the sample with the lid and press MEASUREMENT .

NOTICE

Make sure that the correct sample holder is in the correct orientation for reliable measurement results. For further guidance see Section 3.4.

The result, measurement value and turbidity, will be displayed with the corresponding unit.

6.4 PROCEDURE FOR CORRECT VALUES

To obtain correct values from a measurement, it is crucial to follow the outlined steps:

- Ensure that the cuvette slot is installed correctly (Section 3.4) and that it is kept clean.
- Ensure that the cuvette is free from defects and wiped off in a thorough manner.
- Ensure that the method instructions have been carefully followed.
- Ensure that sample lid is used to cover the sample during measurement, according to instructions.

- Ensure that quality checks are made using calibration standards for both absorbance and turbidity and that method calibration curves are performed regularly.

- Ensure that your IC-VIS device is operating according to specifications and if malfunction is observed, contact Spec-Imaging for technical support.

CAUTION

Spec-Imaging cannot guarantee the safety of the operator if not following the instructions of the chemical method, indicated by the method instruction.

6.5 IC – INSTANT CLARITY TURBIDITY COMPENSATION

Performing colorimetric analysis directly in turbid liquid with a conventional colorimetric system is problematic as measurement results will be biased due to scattering.

By enabling the iC turbidity compensation, the result displayed is compensated reducing the bias caused by turbidity. The iC feature can be enabled/disabled by entering MENU → Settings → Instant Clarity: ON/OFF. The feature status is indicated in the Main Screen according to the icon shown below in Figure 8.



Figure 8. The iC Instant Clarity icon, iC enabled to the left, iC disabled to the right.

NOTICE

Upon delivery, the IC-VIS is pre-configured with a standard iC turbidity compensation setting optimized for the majority of wastewater samples commonly encountered in typical applications. For custom applications and to adapt to other specific water matrices, contact Spec-Imaging support for further instructions.

6.6 TURBIDITY MEASUREMENT

The turbidity measurement is done at any wavelength, in transmission configuration. The device is calibrated at production using Formazine solutions and gives a turbidity reading which agrees with the NTU reading from a standard turbidimeter when measuring a Formazine sample.

NOTICE

Note that turbidity is measured in a transmission configuration and will hence show different turbidity readings compared to a standard nephelometric turbidimeter if measurements are performed on standards and samples different from Formazine.

6.6.1 Turbidity calibration

CAUTION



Chemical exposure: Spec-Imaging cannot guarantee the safety of the operator if not following the instructions in the safety data sheet of the chemicals which may be used for turbidity calibration.

Turbidity can be re-calibrated with Formazine or other secondary turbidity standard at any time. Visit the IC-VIS configurator for instructions.



IC-VIS Configurator

7 MAINTENANCE

WARNING

Multiple hazards. Maintenance should be done according to the safety precautions described in this section.

7.1 INSTRUMENT CLEANING

NOTICE

Do only use a damp cloth with water and regular dish soap to clean the device exterior.

WARNING



Fire hazard: Do not use flammable detergent or solvents for any cleaning of the device.

WARNING



Chemical exposure: Maintain the cleanliness of the cuvette slot in order to maintain reliable measurement results and to prevent operator from chemical exposure.

Keep the device clean and inspect the cleanliness of the cuvette slot with regular intervals to ensure the prerequisites for accurately performed measurements.

The device can be wiped off with a damp cloth and some regular dish soap. Use a lint-free cloth both for the device and the cuvettes when performing measurements.

If device is exposed to splashes of samples, wipe off immediately.

7.2 STORAGE

Store your device with the unit cover on the cuvette slot, and make sure to always remove the cuvette from the cuvette slot after measurements are completed.

7.3 BATTERY REPLACEMENT

⚠ WARNING	
	Fire hazard: Only use LR6 Alkaline AA or HR6 NiMH AA rechargeable batteries exclusively. Refrain from using other substitutes.

⚠ WARNING	
Explosion risk: Do not intermingle new and used batteries and make sure not to use expired batteries. Do not store unit for long periods with batteries installed.	

8 TROUBLESHOOTING

If your IC-VIS is not functioning as expected and you are in need of technical support, please contact support@spec-imaging.com.

9 SERVICE

9.1 REPLACEMENT PARTS AND ACCESSORIES

WARNING

Using other spare parts than the ones by Spec-Imaging is considered as improper use of the device and voids the warranty of the device.

9.2 DISPOSAL

When this product reaches the end of its useful life, it should be disposed of in accordance with local laws and regulations. To find your nearest designated collection point, please contact your local waste disposal authority.

Spec-Imaging takes back old devices and disposes it free of charge.

Please be aware that incorrect disposal of this waste may result in penalties in accordance with local legislation.

9.3 WARRANTY, LIABILITY AND COMPLAINTS

WARRANTY SUMMARY FOR SPEC-IMAGING PRODUCTS

Warranty Terms:

- **Duration:** 1 year from the date of shipment.
- **Eligibility:** Original purchaser only; non-transferable.



- **Coverage:** Includes defects in materials and workmanship under normal use.
- **Remedies:** Spec-Imaging may repair, replace, or refund the defective product, excluding shipping costs.

Warranty Exclusions:

- Damages from improper installation, misuse, or maintenance.
- Accidents, disasters, or external causes.
- Unauthorized modifications or repairs.
- Non-compliance with product instructions.
- Normal wear and tear.

Claim Process:

- Contact Spec-Imaging for return authorization before sending back any product.
- Product inspection is required to confirm the defect and claim validity.

Limitations:

- Repaired or replaced products are covered for the remainder of the original warranty period.
- Purchaser's legal rights may vary by jurisdiction.

Liability and Indemnification:

- **Limitation of Liability:** Company's liability is limited to the product's purchase price. No coverage for indirect, incidental, or consequential damages.

- **No Extended Warranties:** Only covers defects in materials and workmanship; excludes all other warranties, express or implied.
- **Indemnification:** Purchaser to defend and indemnify Spec-Imaging against any claims arising from product misuse or non-compliance with usage instructions.
- **Governing Law:** Laws of Sweden apply.

This summary is part of Spec-Imaging's Terms and Conditions and acceptance of the product constitutes agreement to these terms.

9.4 CONTACT

If you have questions not answered in this manual, then please contact us, see our details below:

Spec-Imaging AB
Ideongatan 3A
Ideon Science Park
223 70 Lund
Sweden

+46 (0) 76 009 81 05

info@spec-imaging.com

www.spec-imaging.com

9.5 VERSION HISTORY

Manual IC-VIS Colorimeter ENG, V1.0,
June 2024

IC-VIS **SAFETY MANUAL** v1.2

Spec-Imaging 2025

Manual IC-VIS Colorimeter ENG, V1.1, January
2025

Manual IC-VIS Colorimeter ENG, V1.2,
February 2025

